




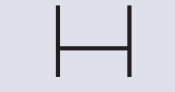

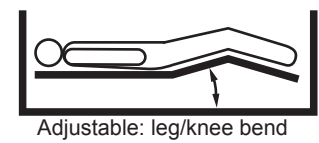
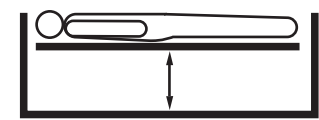
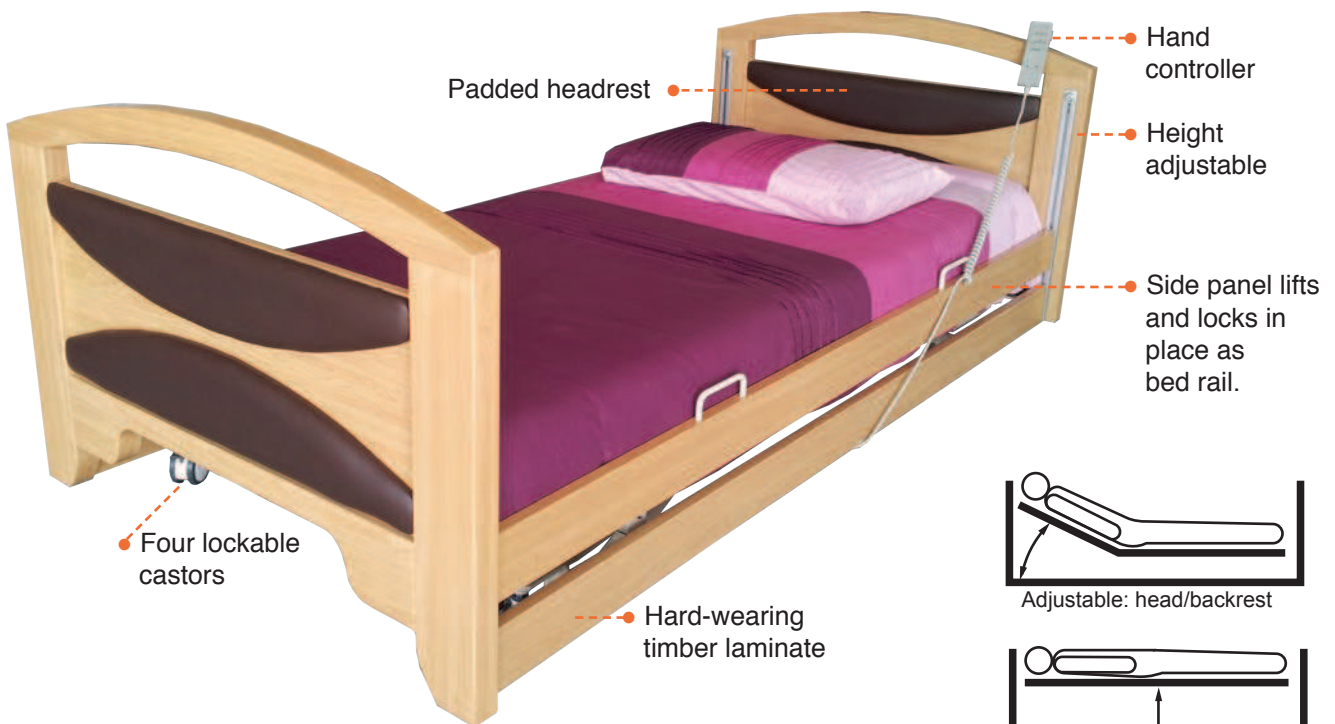


OBM583 MODERNA BED




| Overall Height | Base Height | Overall Length | Overall Width | Castor Size | Unit Weight | Weight Capacity |
|---|--|--|---|---|--|--|
|  74-126cm |  20-72cm |  212cm |  51cm |  7.5cm |  125kg |  180kg |

12 MONTH WARRANTY



OPTIONAL EXTRAS

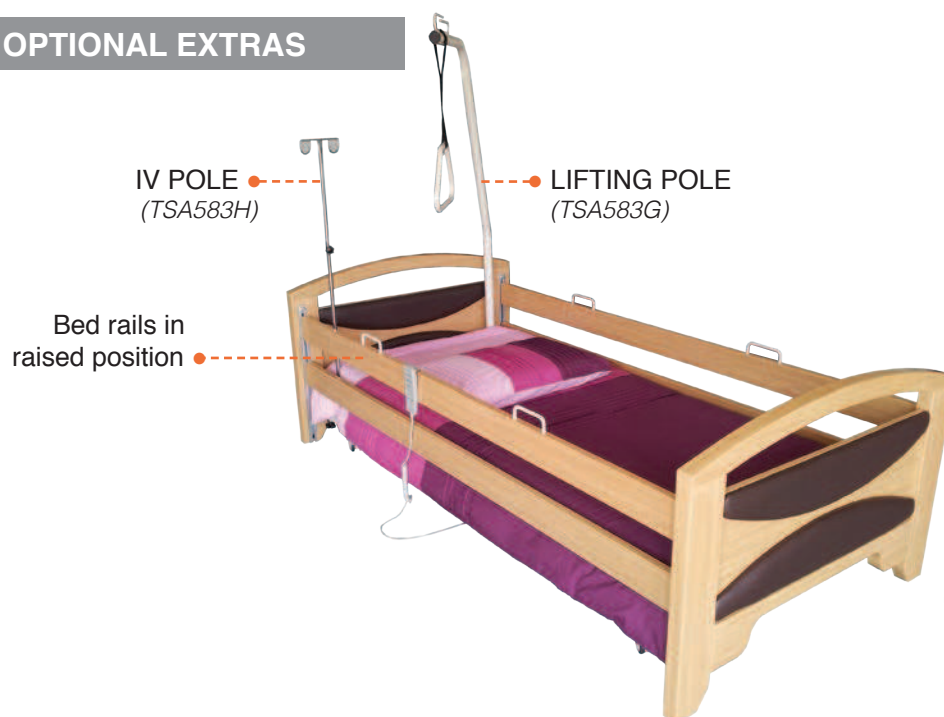

OBM587 MEMORY FOAM MATTRESS HOSPITAL GRADE




| Overall Length | Overall Width | Overall Height | Unit Weight | Weight Capacity |
|----------------|---------------|----------------|-------------|-----------------|
| | | | | |
| 198cm | 90cm | 15.5cm | 25kg | 160kg |

12 MONTH WARRANTY

BEDDING

- Two-way stretch fabric cover
 - removable
 - breathable
 - anti bacterial
 - spill proof
- A top layer of Memory foam contours to the body providing superior comfort and can assist in the prevention of bed sores (bed ulcers).
- A Mattress will contour to the shape of the bed frame when the bed is adjusted into its different positions.

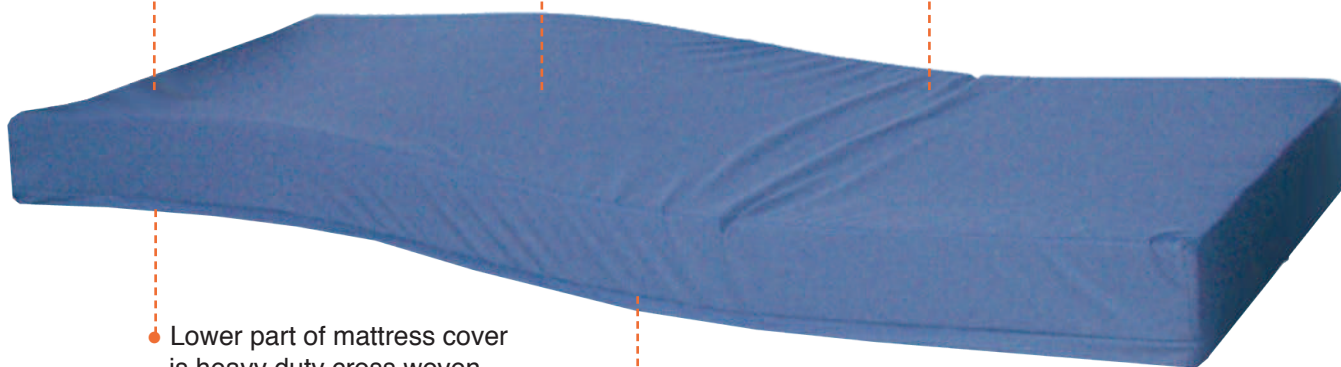

- Lower part of mattress cover is heavy duty cross woven PVC sheet. It will resist tearing if dragged across the floor with the patient lying on top, in the event of an emergency.
- The top and bottom covers are zipped together for easy removal.

- The top layer is memory foam, it has the ability to conform to the shape of the body, providing pressure relief and weight distribution for the patient.



- The sides and bottom of the mattress are made from a firm, high density foam. It provides support and stability to the mattress. The firmer sides create a barrier, reducing the risk of the patient rolling off the mattress.

- The middle layer is a medium density foam, it has a higher resistance to work in combination with the top layer to disperse weight distribution evenly.

Section 7

27

